

Operation Room Tables



Technology for life





Electric C - Arm Compatible Operating Table

USI - 2005

INBUILT KIDNEY BRIDGE

DETACHABLE HEAD SECTION

IMPORTED MOTORS 4 FUNCTION

UNITED **Electric Operation Tables** are designed keeping in view the needs of modern surgery procedures. These tables are versatile & heavy duty enough to cater to all kind of modern surgeries.




CHAIR POSITION



REVERSE ORIENTATION



C - ARM

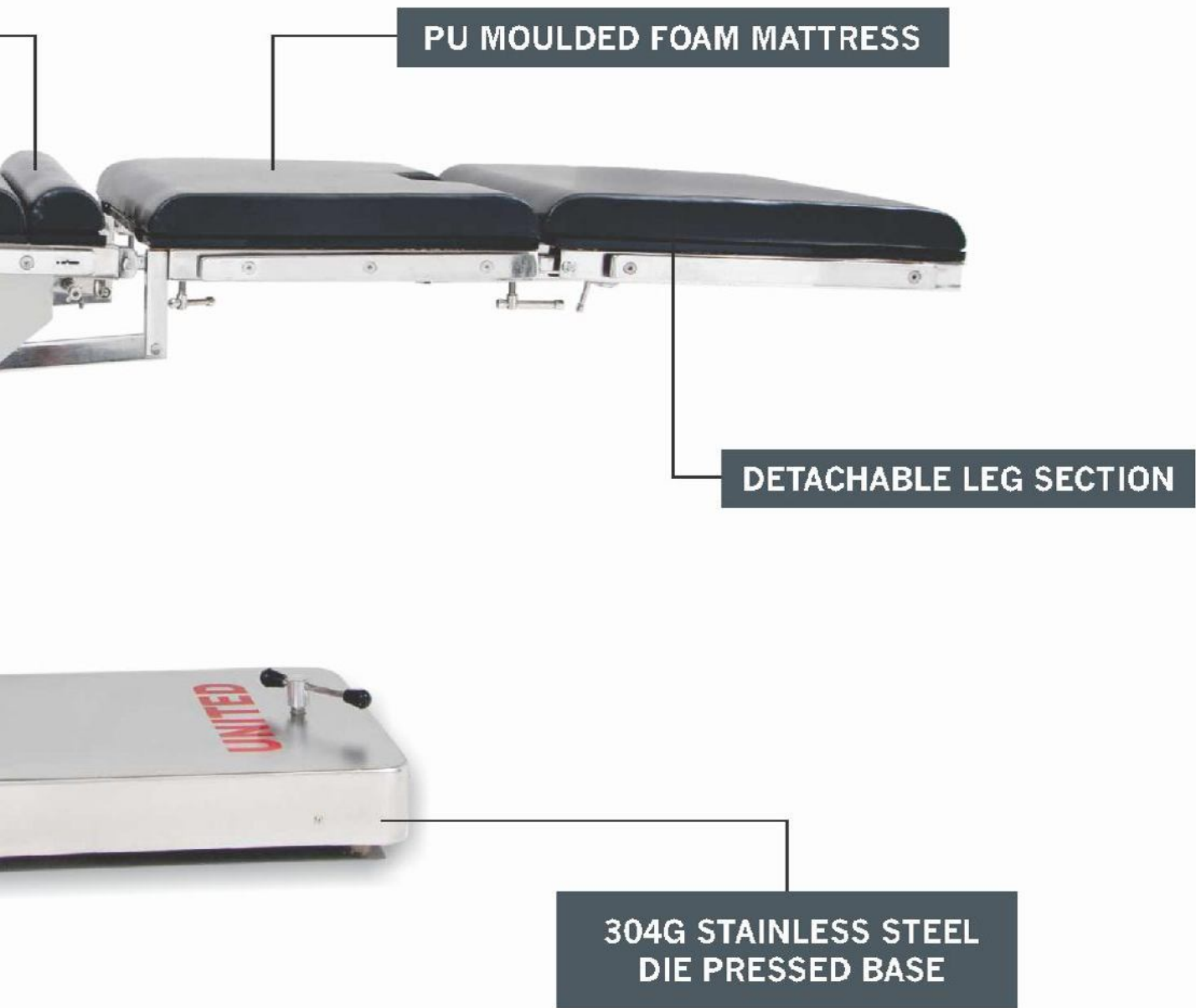

PRONE POSITION



RECTAL SURGERY



LITHOTOMY POSITION

Electric C - Arm Compatible Operating Table

Standard Accessories



Electrical Positions



Tabletop elevation up/down

Mechanical Movements



Leg section adjustment (Up/Down)



Trendelenburg / Ant-Trendelenburg



Flex / Re-Flex



Kidney Bridge



Lateral Tilt Left / Right



Head section adjustment (Up/Down)

Imported Electric linear actuators facilitates the various adjustment of positions. Hi-Lo elevation makes it easier for C-Arm fluoroscopy examination. With built-in kidney bridge suits for kidney operation. Optional traction for orthopedics extend functions of the table. 304G Stainless Steel die pressed cover along with the main body for anti-rust and corrosion.

Hydraulic C - Arm Compatible Operating Table

USI - 2003



Mechanical Movement

1. Flex / Reflex
2. Lateral Tilt
3. Trendlenburg / Reverse Trendlenburg
4. Leg Rest
5. Head Rest
6. Kidney Bridge

Salient Features

- Five section table top allow ease and comfort for surgical positions.
- Radio-translucent top provides easy convenience for X-ray procedures.
- Eccentric Pillar to ensure maximum coverage of C-Arm machine.
- Interchangeable Head & Leg Section.
- Base & column covered with 304G Stainless Steel top to ensure high durability hygiene for the operation room.

Standard Accessories

1. Five sections Mattress :1Set.
2. Anesthetist Screen :1Pcs.
3. Padded Shoulder Support :2Pcs.
4. Padded Lateral Support : 2Pcs.
5. Padded Lithotomy Crutches :2Pcs.
6. Padded Arm Rest :2Pcs.

Electrically Hi-Low

USI - 2003E



Hydraulic Operating Table

Operation Table

Side end Control

USI - 2001



Salient Features

- Five section table top gives ease for various surgical Positions.
- Back section, Trendlenburg/RTB, Flex/Reflex, Kidney section operated through drive handle.
- Hi-Low Position is achieved by concealed hydraulic Pump.
- Base & column covered with 304G Stainless Steel Top to ensure high durability & hygiene for the operation room.
- Leg Section & Head Section Interchangeable.

Standard Accessories

1. Five sections Mattress :1Set.
2. Anesthetist Screen :1Pcs.
3. Padded Shoulder Support :2Pcs.
4. Padded Lateral Support : 2Pcs.
5. Padded Lithotomy Crutches :2Pcs.
6. Padded Arm Rest :2Pcs.
7. Drive Handle :1Pcs.



Specification/Model No.	USI - 2001	USI - 2001E
Overall Length of the table	1950mm	1950mm
Overall width of the table	500mm	500mm
Minimum height of the table	780mm	780mm
Maximum height of the table	1030mm	1030mm
Trendlenburg/Reverse	30degree	30degree
Flex/reflex	80°/220°	80°/220°
Side Tilt	22°	22°
Kidney Section	125mm manually	125mm manually
Leg Plate	90° Down	90° Down
HI-Low lift	Hydraulic	Electric

Orthopedic Attachment

USI - 2006



Orthopedic Attachment

UNITED **Orthopedic Attachment** is a flexible, easily operable instrument for ORTHOPEDIC SURGERY.


The attachment allows the positioning of the extremity as required in regard to the surgery practice and therefore is especially adapted to ensure safe and efficient.



Tibia & Fibula Position

Standard Accessories

- | | |
|--|-----------|
| 1. Orthopedic Section | :1Pc. |
| 2. Pelvic Rest in 2 sizes | :1Pc each |
| 3. Pelvic Support (Radio Translucent) | :1Pc. |
| 4. Inner thigh support (Radio Translucent) | :1Pc. |
| 5. Foot Plate with Ortho shoes | :1Pair |
| 6. Tibia Attachment | :1Pc. |
| 7. Head Surgery table (Radio Translucent) | :1Pc. |
| 8. Steinmen Pin Holding Clamp | :1Pc. |
| 9. Leg holder (Goepel Type) | :1Pair |

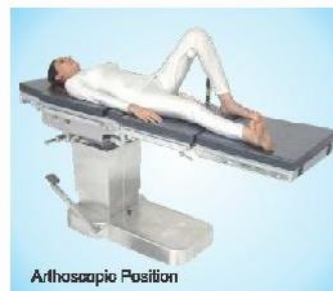

Arthroscopic Position



Neck of Femur Position

Optional

1. Spine Frame Attachment
2. Hand traction Device



Spinal Attachment (Optional)



Hand Surgery

● **Easy Fitting with interlocking facility**

● **304G Stainless Steel**

● **Compatible with USI - 2005, 2003, 2003E**

Technical Data

Specification / Model No.	USI - 2005	USI - 2003	USI - 2003E
Overall Length of the table	1950 mm	1950 mm	1950 mm
Overall width of the table	530 mm	530 mm	530 mm
Minimum height of the table	780 mm	780 mm	780 mm
Maximum height of table	1080 mm	1030 mm	1030 mm
Trendlengburg/Reverse	40°	40°	40°
Flex/Re-flex	80° / 220°	80° / 220° ^c	80° / 220°
Side Tilt	22°	22°	22°
Kidney Section	125 mm manual y	125 mm manual y	125 mm manual y
Leg Section	90° down	90° down	90° down
Head Section	90° down 70° up	90° down 70° up	90° down 70° up
HI-Low Lift	Electric	Hydraulic	Electric



+91 9925063777 | +91 79 6540110
 sales.india@aussin.co
 www.aussin.co | www.ashbme.com



United Surgical Industries

An ISO:9001/13485/14001 & CE Certified Company

Corporate Office:

445, F.I.E. Patparganj Industrial Area, Delhi-110092 (INDIA)
 Ph.: +91-11-22149600 / 94 / 95
 Fax: +91-11-43852679
 E-mail: info@unisurg.com
 Web.: www.unisurg.com

Works:

Pioneer Complex, 12/77-80, Site-IV, Sahibabad Industrial Area
 Distt. - Ghaziabad (U.P.) 201010



Text and image file of this catalogue is property of United Surgical Industries, and protected by copyright and other restriction as well.
 © 2013, United Surgical Industries. All Rights reserved.

The company reserve the right to alter specification and design without our prior notice for the further improvement.